




Disposable Non-Stick Monopolar Suction Coagulator

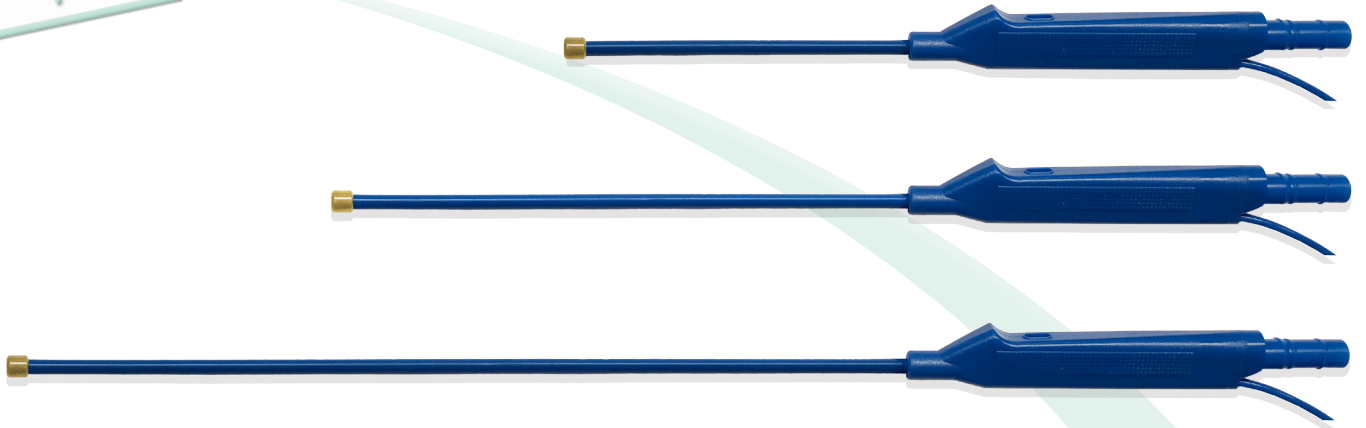

Kirwan Surgical Products disposable non-stick monopolar suction coagulator with resectoscope pin and ferruled tip is offered in three lengths and two French sizes. This non-stick monopolar suction coagulator will allow a surgeon to achieve the same tissue coagulation as a standard monopolar device at lower power settings. This in turn generates less heat, resulting in less thermal damage to the tissue, faster patient recovery time and lowers the risk of post-operative complications.

Kirwan's non-stick monopolar suction coagulator is made with copper tubing, is gold plated and malleable. Similar to the Kirwan AURA technology, the copper tube is significantly more electrically and thermally conductive, allowing surgeons to employ less energy and remove heat faster from tip. In addition, the non-stick monopolar suction coagulator has a bulbous tip, enhancing the non-stick effect and allowing heat to rapidly diffuse from the point of contact, further enhancing the non-stick effect.

30-5202AU Disposable 5 ½" (14cm) Non-stick Monopolar Suction Coagulator with Resectoscope Pin, Ferrule Tip, 10 French

30-5203AU Disposable 8" (20cm) Non-stick Monopolar Suction Coagulator with Resectoscope Pin, Ferrule Tip, 10 French

30-5204AU Disposable 3 ½" (9cm) Non-stick Monopolar Suction Coagulator with Resectoscope Pin, Ferrule Tip, 8 French