



# WE ARE GOING GREEN

Kirwan Surgical's Microdissection Needles are SO popular they deserve their own color!




You will see **GREEN** in your pockets!

Exceptionally priced, Kirwan's distinguished and highly popular microdissection needles are now GREEN in color. By design, Kirwan Surgical's microdissection needles are effective at focusing monopolar energy to the target tissue thereby achieving a precise cut with minimal burning or charring and with minimal damage to peripheral tissue.

The tungsten needle will withstand the rigors of monopolar electrosurgery while maintaining an extremely sharp point. The needle is insulated with a high dielectric coating confining the energy to the tip. The tungsten needle is set into a stainless steel holder which fits a standard cutting/coagulation handpiece.

Go **GREEN**. Go Kirwan **GREEN**.

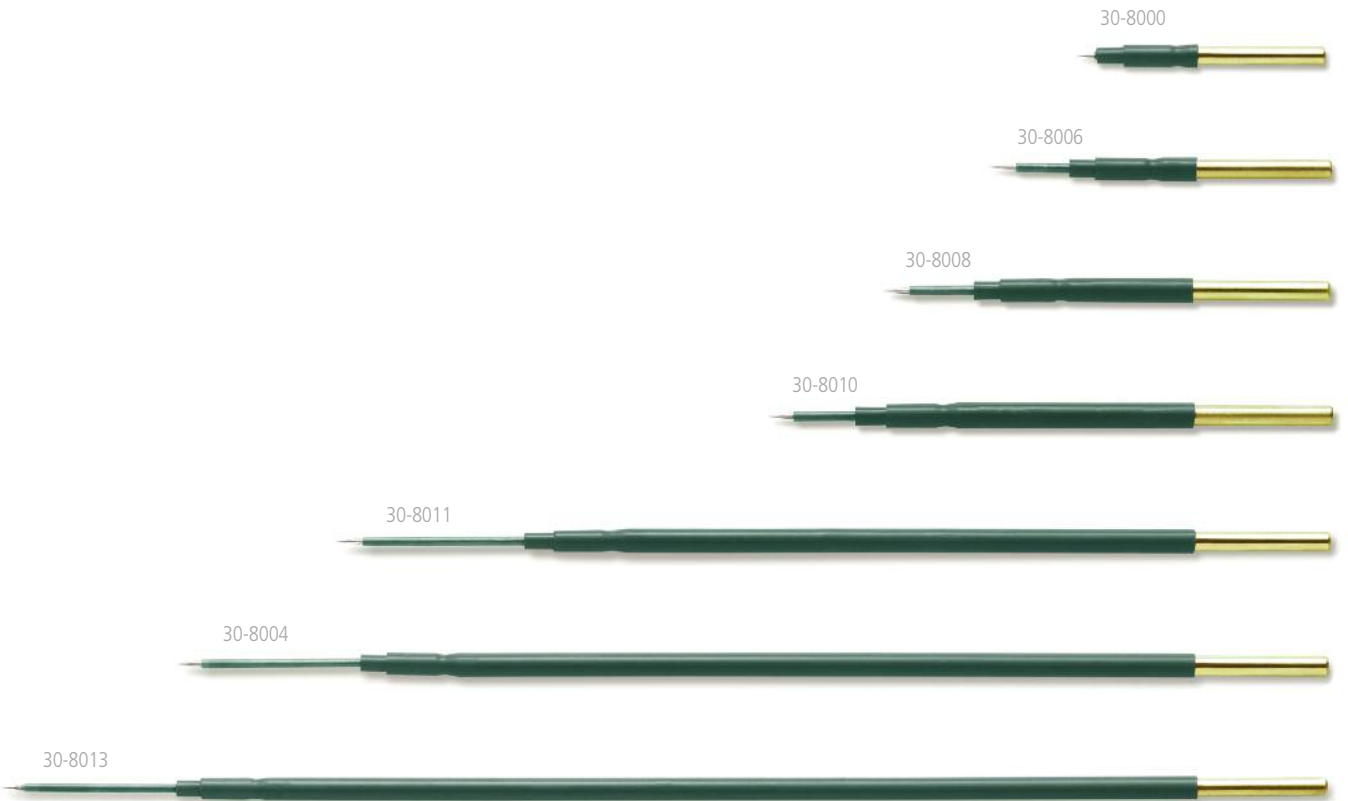

ITEM #	TOTAL LENGTH	DESCRIPTION	FEATURES
30-8000	1 1/4" (3.2cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8006	1 5/8" (4.1cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8008	2" (5.1cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8010	2 1/2" (6.4cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8011	4 5/8" (11.8cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8004	5 1/2" (14cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2
30-8013	6 5/8" (16.8cm)	Disposable Tungsten microdissection needle, insulated, straight	IN 2



ITEM #	TOTAL LENGTH	DESCRIPTION	FEATURES
30-8015	4" (10.2cm)	Disposable Tungsten microdissection curved needle, insulated	IN 2

MICRODISSECTION NEEDLES

# MICRODISSECTION NEEDLES



ITEM #	TOTAL LENGTH	DESCRIPTION	FEATURES
30-8007	1 5/8" (4.1cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN
30-8009	2" (5.1cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN
30-8012	4 1/2" (11.4cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN
30-8014	6 1/2" (16.5cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN
30-8017	4 7/8" (12.4cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN
30-8018	6 7/8" (17.5cm)	Disposable Tungsten microdissection needle, insulated, angled 45 degrees	IN






ITEM #	TOTAL LENGTH	DESCRIPTION	FEATURES
30-8016	5 7/8" (14.9cm)	Disposable Tungsten microdissection curved needle, insulated	IN

MICRODISSECTION NEEDLES



180 Enterprise Drive • Marshfield, MA 02050 USA  
TEL: 781-834-9500 • 888-547-9267 • FAX: 781-834-0022  
WEB: [www.ksp.com](http://www.ksp.com) • EMAIL: [info@ksp.com](mailto:info@ksp.com)

---

Copyright: Kirwan, , WIDESTROKE,  and  are registered trademarks of Kirwan Surgical Products LLC. Products protected under one or more of the following patents, 5,196,009, 5,089,002, 5,290,285, 5,116,333, 5,133,714, 5,989,249, 5,957,928, 4,890,610, 6,749,610, 6,298,550, 6,059,783, 6,228,084, 6,174,310, 6,406,476, 7,621,911, 7,402,754. All Rights Reserved. No part of this catalog may be reproduced or used without written permission. Copyright 2010 Kirwan Surgical Products LLC. All registered trademarks referenced herein are the property of their respective owners.